

LIFEMED T6200 Cardiac

This table features a motorised height adjustment for both patient and practitioner comfort and safety. It also features an electronically operated head section for ease of positioning with the minimum of effort for all practitioners.

The access flap has a simple release mechanism to allow it to drop away 90 degrees for accurate apical imaging across the full range of adult patient sizes. The table top cushioning has been developed and constructed for maximum patient comfort.

The table comes complete with four lockable, rotating castors for ease of movement along with stability and security once in position. The table has optional drop away/removable safety rail.



Specifications

Minimum Height	51cm
Maximum Height	97cm
Cushion Width	70cm
Head Section Length	56cm
Leg Section Length	140cm
Total Length	202cm
Head Tilt Angle	0°/90°
Lift Mechanism	Electric
Breathing Hole	0
Wheels	4 Braked
User Weight Rating	205kg

Warranty

Frame	3 Years
Electronic Components	3 Years
Upholstery	3 Years

Dimensions:


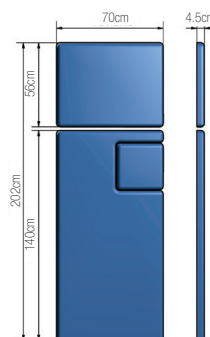


Table Positions:



Maximum
User Weight

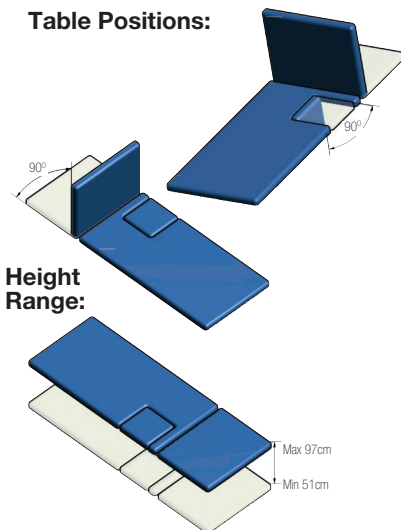

Code	Description
MMC-METCARDIAC	T6200 Cardiac Table - 70cm wide

Optional Accessories:

STOCK	
MMC-METPAPERHOLD	Paper Towel Holder & Towel Roll
NON-STOCK	
MMC-T6006	Cardiac Safety Rail - 89cm