

LIFEMED

T8727 Examination Chair MK1

The Examination Chair MK1, offers the practitioner a variable height chair that can double as a plinth for prone positioning. The chair has swing away arm rests that lift up to allow easy access to the chair. These arm rests are also removable to allow even greater flexibility of use. Each section of the table can be individually moved to a preferred angle to suit treatment requirements.

The backrest is controlled by gas struts giving smooth yet rigid control of this section. A quiet electric actuator is used to raise and lower the chair. A foot switch is provided as standard.

The chair also has detachable arms. Standard with face hole and plug. Comes with four wheels, as standard.



Specifications

Minimum Height	53cm
Maximum Height	95cm
Cushion Width	55cm
Head Section Length	74cm
Mid Section Length	48cm
Leg Section Length	50cm
Total Length	180cm
Total Length (with Tray)	202cm
Head Tilt Angle	0°/80°
Mid Tilt Angle	0°/15°
Leg Tilt Angle	-30°/30°
Lift Mechanism	Electric
Breathing Hole	1
Retractable Wheels	Standard x 4
User Weight Rating	205kg

Warranty

Frame	3 Years
Electronic Components	3 Years
Upholstery	3 Years

Dimensions:



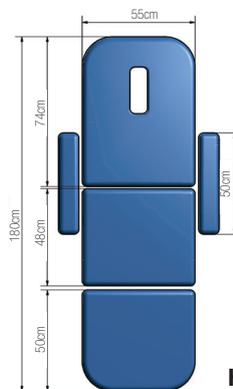


Table Positions:

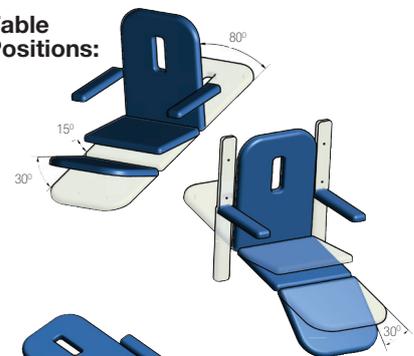

Height Range:

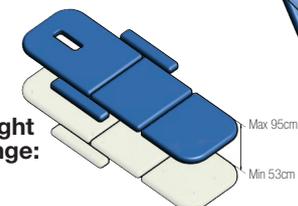

Code	Description
MMC-METPOD	T8727 Examination Chair MK1 - 55cm wide